


TIR-1™

Non-contact, clinical-grade infrared thermometer with Bluetooth® connectivity




Non-contact module reduces the risk for patient cross-contamination while reducing costs and waste by eliminating the need for probe covers and other disposables

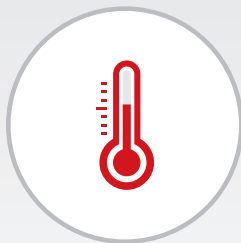

Clinical-grade, infrared thermometer provides forehead temperature measurement across all patient populations



Simple, one-button operation delivers results in seconds, reducing measurement time compared to traditional methods

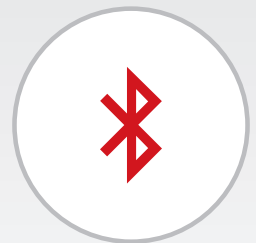

Bluetooth automates data transfer to a connected Masimo device, enabling streamlined integration into the bedside device and EMR

Simplifying Workflows with Automated Electronic Charting



Temperature measurement is transferred to the Masimo device via Bluetooth



Patient data, including temperature, are aggregated by the Masimo device



Masimo Patient SafetyNet™* or Iris Gateway™ converts data into HL7



Patient SafetyNet or Iris Gateway automates data transfer to the hospital EMR

Integration

TIR-1 integrates with the following Masimo devices:



Root®



Root with Noninvasive Blood Pressure and Temperature

TIR-1 Specifications

| | |
|-----------------------------------------------------------------------------|----------------------------------------------------------------------|
| ACCURACY | ENVIRONMENTAL |
| Display Range 34.4° - 42.2° C (94° - 108° F) | Storage Temperature -20° - 60° C (-4° - 140° F) |
| PHYSICAL CHARACTERISTICS | Operating Humidity 15 - 85% relative humidity (non-condensing) |
| Power 2x 1.5V, AA Alkaline Batteries | COMPLIANCE |
| Dimensions 150 mm x 48.48 mm x 55.16 mm (5.90 in x 1.91 in x 2.17 in) | ASTM E1965-1998 |
| ORDERING INFORMATION | This product meets the requirements of the applicable EC directives. |
| Part Number 4301 | |

* The use of the trademark SafetyNet is under license from University HealthSystem Consortium

Root with TIR-1 is not licensed for sale in Canada.

For professional use. See instructions for use for full prescribing information, including indications, contraindications, warnings, and precautions.

Masimo U.S.
Tel: 1 877 4 Masimo
info-america@masimo.com

Masimo International
Tel: +41 32 720 1111
info-international@masimo.com

