

Halo ION™

A Comprehensive, Personalised, Continuous Patient Score

Halo ION allows clinicians to aggregate trend data from multiple physiologic parameters into a single, continuous patient score which is configured according to hospital protocol.




Halo ION

- > Creates a continuous patient score according to hospital protocol by examining key characteristics of multiple parameters for more personalised care
- > Displays the relative contribution of each parameter, helping clinicians better assess possible causes of patient deterioration as defined by hospital protocol
- > Automates the process by which clinicians assess patient status, reducing cognitive overload and streamlining workflows


Automate Workflows with Halo ION

Halo ION integrates values from hospital protocol, including acceptable parameter limits and relative parameter importance, to automatically calculate a continuous, cumulative score, which may help clinicians assess patient status.

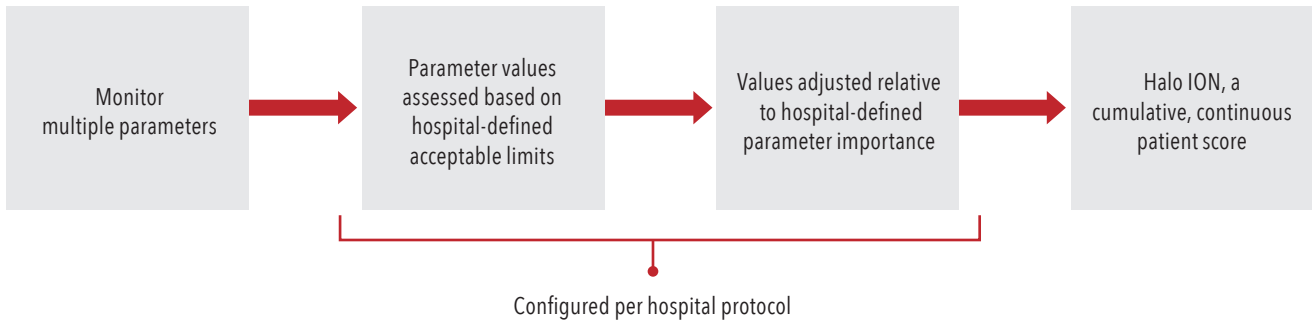

A Comprehensive, Personalised, Continuous Patient Score

The relative contribution of each parameter to the Halo ION score is shown by the size of the bubble next to each parameter, which may help clinicians assess possible causes of patient deterioration per hospital protocol

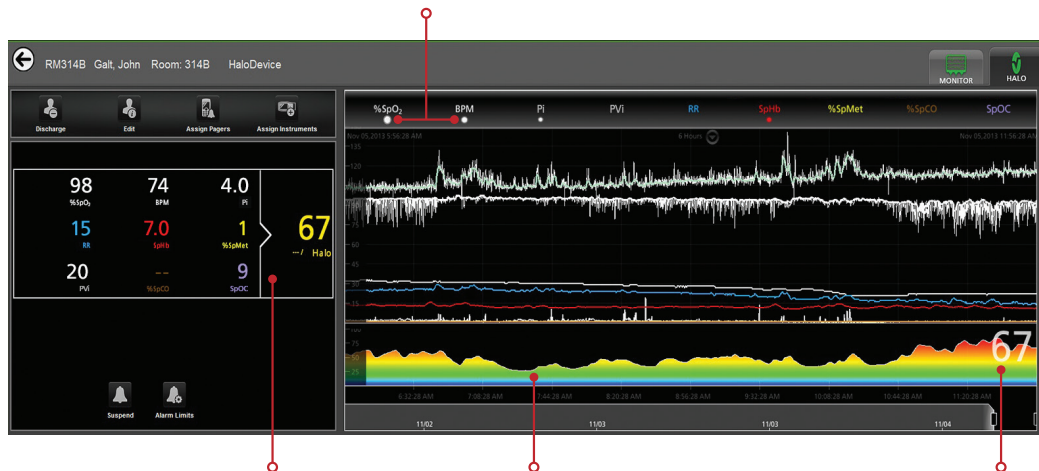

Aggregates data from multiple physiologic parameters into one score per hospital protocol

Continuous trend data displays relative changes in Halo ION score over time

Displays a single score ranging from 0 to 100; a high score is indicative of patient risk per hospital protocol

Specifications

Patient SafetyNet™* 9095 Patient SafetyNet Service 119-XY

* The use of the trademark PATIENT SAFETYNET is under license from University HealthSystem Consortium.

For professional use. See instructions for use for full prescribing information, including indications, contraindications, warnings, and precautions.

Masimo U.S.
Tel: 1 877 4 Masimo
info-america@masimo.com

Masimo International
Tel: +41 32 720 1111
info-international@masimo.com

