


Early Warning Score

Automatic Device-based Calculation



Spot-check Calculation with Time/Date Stamp

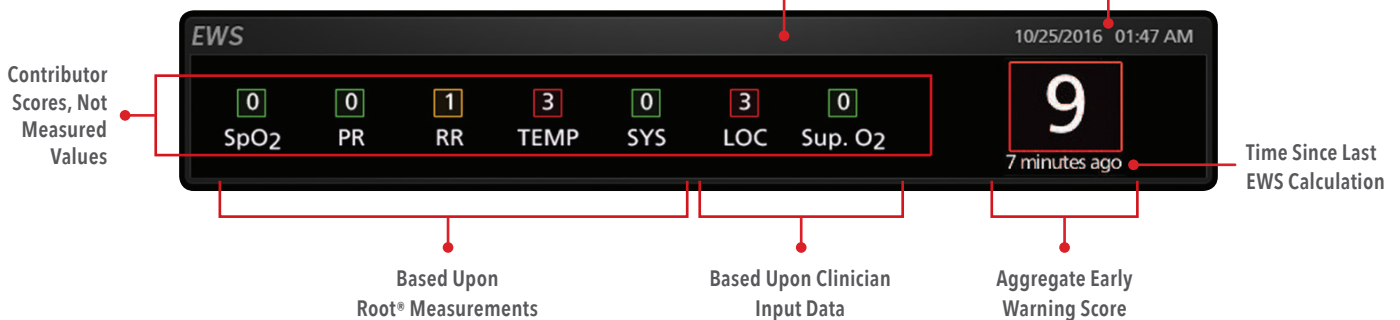


- > Early Warning Score (EWS) is initiated by the clinician and automatically calculated
- > EWS contributor scores are calculated using measured values and clinician input, then combined into an aggregate EWS
- > Root features an instantly interpretable, high-visibility display with intuitive, multi-touch navigation for easy and adaptable use in hospital environments



Advanced Connectivity Solutions with Patient SafetyNet™* or Iris Gateway™

Bedside Patient Association

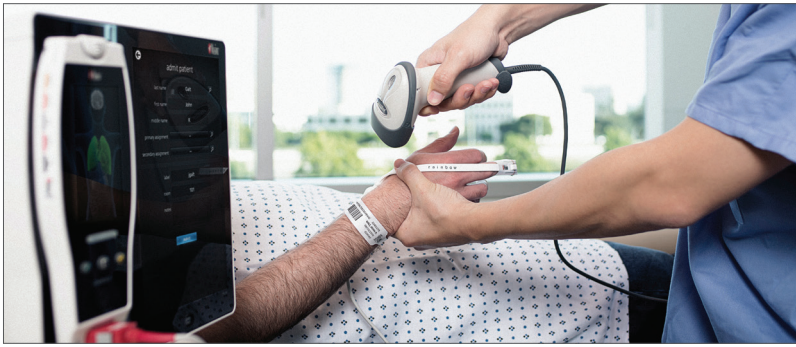

Interface with hospital HL7 Admit, Discharge, and Transfer (ADT) system to associate the patient, device, and clinician at the bedside

Automated Electronic Charting



Patient data, including EWS, is captured by Root




Patient SafetyNet or Iris Gateway converts data into HL7



Patient SafetyNet or Iris Gateway automates data transfer to the hospital EMR

Manual Electronic Charting at the Bedside

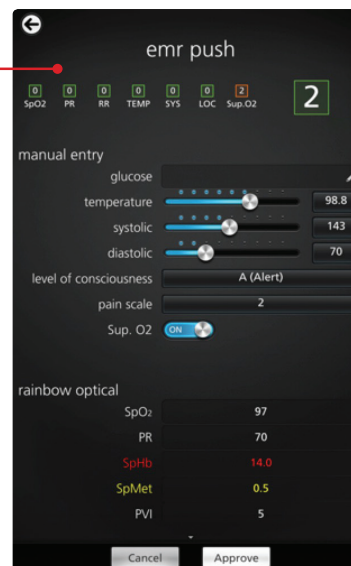

To simplify documentation of patient data, clinicians can also manually enter and verify data on Root, including level of consciousness, pain scale, and urine output, and send to the EMR

* The use of the trademark Patient SafetyNet is under license from University Health System Consortium

For professional use. See instructions for use for full prescribing information, including indications, contraindications, warnings, and precautions.

Masimo U.S.
Tel: 1 877 4 Masimo
info-america@masimo.com

Masimo International
Tel: +41 32 720 1111
info-international@masimo.com

